BIOS living




18”/45 cm Wheelchair Instruction Manual 56084

Thank you for purchasing the BIOS Living wheelchair. Please review the manual before using. **DO NOT assemble or use this equipment without first reading and understanding this instruction manual.**

Note: This owner’s manual was compiled from the most recent specifications and product information available at the time of publication. Thermor Ltd. reserves the right to make changes, which may cause slight variations between the illustration and explanations in this manual and the product purchased.

FEATURES

- Seat width: 18”/46 cm
- Folds easily
- 8”/20 cm front wheels
- 24”/61 cm rear MAG wheels
- Dual locking rear wheels for safe and secure transfers
- Adjustable, detachable footrests
- Comfortable handle grips offer increased comfort and control
- Convenient storage pocket
- Fully padded seat and backrest with all-weather easy to clean nylon covering

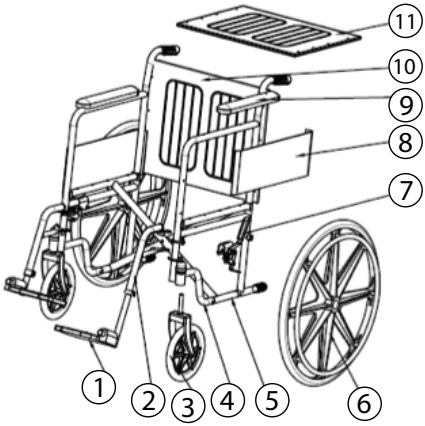

1. Footrest
2. Metal foot tube
3. 8”/20 cm front wheels
4. Cross support
5. Side support bars
6. 24”/61 cm MAG rear wheels
7. Brake
8. Metal side plate
9. Armrest
10. Backrest
11. Padded seat

Assembly Instructions

1. Unfold the wheelchair

Push both sides of the seat down to unfold the wheelchair to normal operating condition.

2. Assemble the footrests

Assemble both footrests to the front end of the wheelchair as shown in the diagram. Move and turn both left and right footrests to the normal operating position (facing forward).


3. Fold the wheelchair

Turn the footrests to the vertical position and pull-up on the middle of the seat to fold the wheelchair.

⚠ Safety Precaution



- The maximum load for the wheelchair is 300 lbs/136 kg, this includes the rider and any articles he/she may be carrying.
- Drive the wheelchair on smooth ground, and avoid obstacles. Failure to do so may cause the wheelchair to overturn.
- **Do not** use the wheelchair on stairs or escalators.
- **Do not** use the wheelchair as a seat in a motor vehicle, unless the motor vehicle is equipped with certified equipment to hold and support the wheelchair.
- **Do not** stand on the footrests, standing on the footrest could cause the wheelchair to overturn

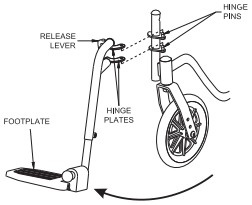



Maintenance

- Check if there are loose bolts before operation and/or within the first month of operation. Tighten loose bolts immediately.
- Check every three months during normal operation, to ensure that all parts are in good condition.
- Check all nuts of the wheelchair (especially the fixing nuts on the rear wheel axle). Adjust promptly, if loose. The connection bolts on the wheelchair cross bars are designed to be loose - no adjustments are necessary.
- Regularly check the movable parts and turning mechanism of the wheelchair, and apply lubricant.
- Dry the wheelchair promptly if it gets wet.
- Regularly wipe the wheelchair with a soft cloth during normal operation and apply antirust wax.

BIOS Living Warranty for Wheelchairs

100% Satisfaction Guarantee

If at any time, you are not completely satisfied with the performance of this device, call our BIOS Medical Hotline and speak with a customer service representative, who will make arrangements to have the device tested or replaced to your full satisfaction.

If you have questions regarding the operation of your wheelchair call the **BIOS Medical customer service: 1-800-387-8520.**

Warranty Coverage 3 Years: Frames are warrantied for 3 years against defects in materials and/or workmanship for the following models: Phoenix Series wheelchairs 56184, 56185, 56186 and BIOS Series Wheelchairs 56084

Warranty Coverage 1 Year: One year warranty on the parts and upholstery to be free from defect in materials and workmanship, excluding normal wear and tear.

This warranty does not cover device failure due to owner misuse or negligence, or normal wear and tear.

This warranty extends only to the original purchaser and is not transferrable. Thermor Ltd. will not warranty any parts that have been damaged due to improper installation, modification and/or misuse.

It is the sole discretion of Thermor Ltd. to determine best method of repair/replacement. The replaced part will only be replaced once under the original warranty.

This warranty does not cover any transport/shipping costs. This warranty does not apply if the product is subject to misuse, neglect, rough handling or damage. Ship the unit prepaid and insured (at owners option) to:

Thermor Ltd.
Returns Department
16975 Leslie Street
Newmarket, ON L3Y 9A1
www.biosmedical.com
Email: support@biosmedical.com

Please include your name, return address, phone number, and email address. Thermor will test or replace (at Thermor's option) free of charge any parts necessary to correct the defect in material or workmanship.

Please allow 10 days for return shipping.

Paper
Papier
BIOS | Medical
MANUFACTURED BY/FABRIQUÉ PAR :
THERMOR LTD.
16975 LESLIE STREET
NEWMARKET, ON L3Y 9A1
MADE IN CHINA/FABRIQUÉ EN CHINE
WWW.BIOSMEDICAL.COM
100% Canadian Owned & Operated
Propriétaire-exploitant canadien à 100 %