LIMITED WARRANTY

The Fingertip Pulse Oximeter has a 2 year warranty to be free of manufacturing defects in materials and workmanship under normal applications for the original owner. If this product becomes inoperable due to defect and requires correction, return the product with all component pieces and proof of purchase to the address listed below. This warranty does not cover any shipping/transport costs. This warranty does not apply if the product is subject to misuse, neglect, rough handling or damage.

Ship the unit prepaid and insured (at owner's option) to:
Thermor Ltd.
Attn: Returns Department
16975 Leslie Street
Newmarket, ON L3Y 9A1

www.biosmedical.com
Email: support@biosmedical.com

Limitations

This warranty covers all defects encountered in normal use of the equipment and does not apply in the following cases.

1) If the equipment has been serviced by other than a certified center.

2) Damage to the equipment due to mishandling, abuse, accident or not following operating instructions.

3) INSTRUMENTS: Warranty does not extend to display face, batteries, carrying case or lanyard.

For troubleshooting assistance please contact Consumer Support. Our representatives are trained to provide you assistance over the phone and may be able to resolve your problem without returning. For the best possible assistance please have your unit available when calling us.

Service Phone: 1 (866) 517-2970

www.biosmedical.com

OPERATOR'S MANUAL

MD300C29

Fingertip Pulse Oximeter

General Description

Oxygen binds to hemoglobin in red blood cells when moving through the lungs. It is transported throughout the body as arterial blood. A pulse oximeter uses two frequencies of light (red and infrared) to determine the percentage (%) of hemoglobin in the blood that is saturated with oxygen. The percentage is called blood oxygen saturation, or SpO_2. A pulse oximeter also measures and displays the pulse rate at the same time it measures the SpO_2 level.
Diagram of Operation Principle
1. Red and Infrared-ray Emission Tube
2. Red and Infrared-ray Receipt Tube

Precautions For Use

1. Please read the manual carefully before use.
2. Operation of the pulse oximeter may be affected by the use of an electrosurgical unit (ESU).
3. The pulse oximeter must be able to measure the pulse properly to obtain an accurate SpO_2 measurement. Verify that nothing is hindering the pulse measurement before relying on the SpO_2 measurement.
4. Do not use the pulse oximeter in an MRI or CT environment.
5. Do not use the pulse oximeter in situations where alarms are required. The device has no alarms. It is not for continuous monitoring.
6. Do not use the pulse oximeter in an explosive atmosphere.
7. The pulse oximeter is intended only as an adjunct in patient assessment. It must be used in conjunction with other methods of assessing clinical signs and symptoms.
8. In order to ensure correct sensor alignment and skin integrity, the maximum application time at a single site for our device should be less than half an hour.
9. Do not sterilize the device using autoclaving, ethylene oxide sterilizing, or immersing the device in liquid. The device is not intended for sterilization.
10. Follow local ordinances and recycling instructions regarding disposal or recycling of the device and device components, including batteries.
11. This equipment complies with IEC 60601-1-2 for electromagnetic compatibility for medical electrical equipment and/or systems. However, because of the proliferation of radio-frequency transmitting equipment and other sources of electrical noise in healthcare and other environments, it is possible that high levels of such interference due to close proximity or strength of a source might disrupt the performance of this device.
12. Portable and mobile RF communications equipment can affect medical electrical equipment. The portable and mobile RF communications equipment should be used no closer than 30cm (12 inches) to any part of the device, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.
13. This equipment is not intended for use during patient transport outside the healthcare facility.
14. The patient is an intended operator. The patient can safely use all functions of the device.
15. It may be unsafe to:
 —use accessories, detachable parts and materials not described in the instructions for use.
 —interconnect this equipment with other equipment not described in the instructions for use.
 —disassemble, repair or modify the equipment.
16. The material that contact with the patient's skin has passed the ISO10993-5 Tests for invitro cytotoxicity and ISO10993-10 Tests for irritation and delayed-type hypersensitivity.
17. Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally.
18. The use of accessories, transducers and cables other than those specified or provided by the manufacturer of this equipment could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation.
19. When the signal is not stable, the reading may be inaccurate. Please do not refer to it.
20. The material of the device has no nature latex.
21. The pulse oximeter equipment is calibrated to display functional oxygen saturation.
22. The waveform we provide is normalized.
23. Do not modify this equipment without authorization of the manufacturer.
24. The manufacturer will make available on request circuit diagrams, component part lists, descriptions, calibration instructions, or other information that will assist the user's appropriately trained personnel to repair those parts of the equipment designated by the manufacturer to be repairable.
25. Stop using and contact local service center if one of the following cases occurs:
 -Any of the problems in the Possible Problems and solutions cannot be solved.
 -The oximeter cannot be powered on in any case and not the reasons of battery.
 -There is a crack on the oximeter or damage on the display resulting readings cannot be identified; the spring is invalid; or the key is unresponsive or unavailable.
26. Aging infrared-ray detector or insufficient battery level may affect the equipment performance. Please follow the instructions in the manual to maintain the device.
27. Please contact the manufacturer for any question about the usage or maintenance.

Inaccurate measurements may be caused by

1.Significant levels of dysfunctional hemoglobin (such as carbonyl - hemoglobin or methemoglobin).
2. Intravascular dyes such as indocyanine green or methylene blue.
3. High ambient light. Shield the sensor area if necessary.
4. Excessive patient movement.
5. High-frequency electrosurgical interference and defibrillators.
6. Venous pulsations.
7. Placement of a sensor on an extremity with a blood pressure cuff, arterial catheter, or intravascular line.
8. The patient has hypotension, severe vasoconstriction, severe anemia, or hypothermia.
9. The patient is in cardiac arrest or is in shock.
10. Fingernail polish or false fingernails.
11. Weak pulse quality (low perfusion).
12. Low hemoglobin.

Product Features

1. High brightness LCD display SpO_2, PR, Pulse bar and Plethysmograph.
2. Low power consumption; battery-low indicator
3. Automatically power off; 2 AAA alkaline batteries
4. When no signal is detected, the pulse oximeter will power off automatically in 8 seconds.

Contraindication

Not found yet.

Intended Use


The Fingertip Pulse Oximeter is a handheld non-invasive device intended for spot-checking of oxygen saturation of arterial hemoglobin (SpO_2) and Pulse Rate of adult, adolescent and child patients in hospitals, hospital-type facilities and homecare.

Operation Instructions

1. Install two AAA batteries into battery compartment correctly.
2. Place clamp over finger nail as the right diagram.
3. Insert one finger into rubber hole of the oximeter fully.
4. Press the switch once on the front panel.
5. Finger and body should not tremble during measuring.
6. Read correspondent data from display screen.
7. Six display modes.

After turning on the oximeter, each time you press the power switch, the oximeter will switch to another display mode. There are 6 display modes shown as follows:

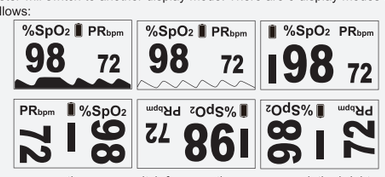When you press the power switch for more than one second, the brightness of the oximeter will be changed by degrees, indicated Br 1 to Br 10. There are 10 levels on brightness; the default is level four.


NOTE: Please use medical alcohol to clean the rubber inside the oximeter, and clean the test finger using alcohol before and after each test. (The rubber inside of the oximeter belongs medical rubber, which has no toxin and no harmful to the skin of human being).

Product Accessories

1. One lanyard
2. Two AAA batteries
3. One user's manual
4. One carrying case

Battery Installation

1. Install two AAA batteries into the battery compartment. Match the plus (+) and minus (-) signs in the compartment. If the polarities are not matched, damage may be caused to the oximeter.
2. Slide the battery door cover horizontally along the arrow .

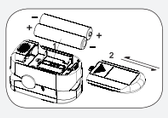Notes:
Please remove the batteries if the pulse oximeter will not be used for long periods of time.

Maintenance and Storage

1. Replace the batteries in a timely manner when low voltage lamp is lighted.
2. Clean surface of the fingertip oximeter before it is used in diagnosis for patients.
3. Remove the batteries if the oximeter is not operated for a long time.
4. It is best to store the product in -25℃ ~+70°C and ≤93% humidity.
5. Keep in a dry place. Extreme moisture may affect oximeter lifetime and may cause damage.
6. Dispose of battery properly; follow any applicable local battery disposal laws.

Cleaning the fingertip pulse oximeter
Please use medical alcohol to clean the silicone touching the finger inside of oximeter with a soft cloth dampened with 70% isopropyl alcohol. Also clean the being tested finger using alcohol before and after each test.
Do not pour or spray liquids onto the oximeter, and do not allow any liquid to enter any openings in the device. Allow the oximeter to dry thoroughly before reuse.
The fingertip pulse oximeter requires no routine calibration or maintenance other than replacement of batteries.
The use life of the device is five years when it is used for 15 measurements every day and 10 minutes per one measurement. Stop using and contact local service center if one of the following cases occurs:
● An error in the Possible Problems and solutions is displayed on screen.
● The oximeter cannot be powered on in any case and not the reasons of battery.
● There is a crack on the oximeter or damage on the display resulting readings cannot be identified; the spring is invalid; or the key is unresponsive or unavailable.


Disinfecting
The applied parts touching the patients' body are required to be disinfected once after each use. The recommended disinfectants include: ethanol 70%, isopropanol 70%, glutaraldehyde-type 2% liquid disinfectants.
Disinfection may cause damage to the equipment and is therefore not recommended for this pulse oximeter unless otherwise indicated in your hospital's servicing schedule. Clean the pulse oximeter before disinfecting it.
CAUTION: Never use EtO or formaldehyde for disinfection.

Specifications

1. Display Type:
Colored display
2. SpO_2
Measurement range: 70%-100%
Accuracy: 70%~100%: ±2%; 0%~69% no definition
Resolution: 1%
Note: A functional tester cannot be used to assess the accuracy of a pulse oximeter monitor or sensor. Clinical testing is used to establish the SpO2 accuracy. The measured arterial hemoglobin saturation value (SpO2) of the sensors is compared to arterial hemoglobin oxygen (SaO2) value, determined from blood samples with a laboratory CO-oximeter. The accuracy of the sensors in comparison to the CO-oximeter samples measured over the SpO2 range of 70%~100%. Accuracy data is calculated using the root-mean-squared (Arms value) for all subjects, per ISO 80601-2-61, Medical Electrical Equipment–Particular requirements for the basic safety and essential performance of pulse oximeter equipment for medical use.
A functional tester is used to measure how accurately Pulse Oximeter is reproducing the specified calibration curve and the PR accuracy.
The model of functional tester is Index2 FLUKE simulator and the version is 2.1.3.
3. Pulse Rate
Measure range: 30bpm~250bpm
Accuracy: 30bpm~99bpm, ±2bpm; 100bpm~250bpm, ±2%
Resolution: 1bpm
4. Probe LED Specifications

	Wavelength	Radiant Power
RED	660±3nm	3.2mW
IR	905±10nm	2.4mW

NOTE: The information about wavelength range can be especially useful to clinicians.
5. Power Requirements
Two AAA alkaline Batteries
Power consumption: Less than 40mA
Battery Life: Two AAA 1.5V, 1200mAh alkaline batteries could be continuously operated as long as 18 hours.
6. Environment Requirements
Operation Temperature: 5°C ~ 40°C
Storage Temperature: -25°C ~ +70°C
Ambient Humidity: 15%~ 93% no condensation in operation;
≤93% no condensation in storage/transport
Atmosphere pressure: 70kPa~106kPa
Notes: When the ambient temperature is 20°C, it is required 6 hours for the equipment to warm from the minimum storage temperature or 4 hours to cool from the maximum storage temperature between uses until it is ready for its intended use.
7. Equipment Response Time
As shown in the following figure.
Response time of slower average is 12.4 s.

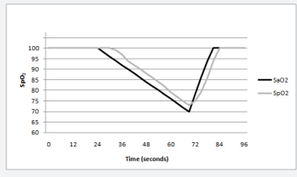
8. Classification
According to the type of protection against electric shock: INTERNALLY POWERED EQUIPMENT;
According to the degree of protection against electric shock: TYPE BF APPLIED PART (applied part: the rubber hole of the device);
According to the degree of protection against ingress of water: IP22;
According to the mode of operation: CONTINUOUS OPERATION.

Clinical Study Summary

The following details are provided to disclose actual performance observed in the clinical validation study of healthy adult volunteers. The ARMS value analysis statement and Bland-Altman plot of data is shown as following:
ARMS Value Analysis Statement

Item	70--100	90--100	80--<90	70--<80	70-<90
#pts	227	91	81	55	136
Bias	0.99	0.73	1.44	0.74	1.16
ARMS	1.80	1.16	2.04	2.21	2.11

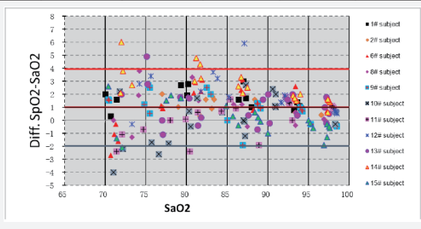Bland-Altman Plot Graphic

Possible Problems and Solutions

Problems	Possible reason	Solution
SpO_2 or PR can not be shown normally	1. Finger is not inserted correctly. 2. User's oxyhemoglobin value is too low to be measured. 3. Cold hands, poor circulation, weak pulse and movement may cause unstable measurement results.	1. Retry by inserting the finger. 2. There is excessive illumination. 3. Please keep your hands still and warm and retry. Try some more times.
SpO_2 or PR is shown unstably	1. Finger might not be inserted deep enough. 2. Finger is trembling or user's body is in movement status.	1. Retry by inserting the finger. 2. Try not to move.
The oximeter can not be powered on	1. Batteries might be installed incorrectly. 2. The batteries' power might be inadequate or not be there at all. 3. The oximeter might be damaged.	1. Please reinstall the batteries according to the polarity marked inside the battery compartment. 2. Please replace batteries. 3. Please contact with local customer service center.
Indication lamps are suddenly off	1. The product is automatically powered off when no signal is detected longer than 8 seconds. 2. The battery power is too low to work.	1. Normal. 2. Replace the batteries.
"Err7" is displayed on screen	Err 7 means all the emission LED or reception diode in damaged.	Please contact with local customer service center.

Symbol Definitions

Symbol	Definition	Symbol	Definition
	Type BF applied part.	SpO_2	No SpO_2 Alarm
	Follow instruction for use		Power switch
IP22	Protected against dripping water	LOT	Lot No.
SpO_2	Oxygen saturation	-25℃ min, +70℃ max, RH≤93% non-condensing	Storage temperature and relative humidity
PR_{bpm}	Pulse Rate (bpm)		Manufacturer's information
	Low power indication		Full Power Indicator
	Attention		Conformity to WEEE Directive
	Date of Manufacture		

Notes:
1. The illustrations used in this manual may differ slightly from the appearance of the actual product.
2. The specifications are subject to change without prior notice.

Electromagnetic Compatibility

The device conforms to IEC60601-1-2 Electromagnetic Compatibility (EMC) standard. Essential performance is defined as SpO2 accuracy and pulse rate accuracy or an indication of abnormal operation. Accuracies may be affected as a result of exposure to electromagnetic disturbances that are outside of the environments listed in the intended use. If issues are experienced, move the device away from the source of electromagnetic disturbances.

Table 1: Electromagnetic Emissions Limits and Compliance

Emission test	Compliance
RF emissions CISPR 11	Group 1
Note: Harmonic Emissions (IEC 61000-3-2), Voltage Flicker Emissions (IEC 61000-3-3) are not applicable.	

Table 2: Electromagnetic Immunity

Immunity Test	Compliance	
Electrostatic Discharge (ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	
Rated power Frequency Magnetic Fields IEC 61000-4-8	30 A/m 50Hz and 60 Hz	
Radiated RF	80 MHz – 2.7 GHz	10 V/m 80% AM 1kHz
	380 – 390 MHz	27 V/m Pulse mod. 18Hz
	430 – 470 MHz	28 V/m FM±5Hz deviation 1kHz sine
	704 – 787 MHz	9 V/m Pulse mod. 217Hz
	800 – 960 MHz	28 V/m Pulse mod. 18Hz
	1.7 – 1.99 GHz	28 V/m Pulse mod. 217Hz
	2.4 – 2.57 GHz	28 V/m Pulse mod. 217Hz
	5.1 – 5.8 GHz	9 V/m Pulse mod. 217Hz
Note: Electrical Fast Transients (IEC 61000-4-4), Surge (IEC 61000-4-5), Voltage dips (IEC 61000-4-11), Conducted Immunity (IEC 61000-4-6) are not applicable.		

Test frequency	Modulation	IMMUNITY TEST LEVEL (A/m)
30 kHz [a]	CW	8
134.2 kHz	Pulse modulation [b] 2.1kHz	65 [c]
13.56 MHz	Pulse modulation [b] 50kHz	7.5 [c]

a) This test is applicable only to ME EQUIPMENT and ME SYSTEMS intended for use in the HOME HEALTHCARE ENVIRONMENT.

b) The carrier shall be modulated using a 50% duty cycle square wave signal.

c) r.m.s., before modulation is applied.

Distributed by:
Thermor Ltd.
16975 Leslie Street
Newmarket, ON L3Y 9A1
WWW.BIOSMEDICAL.COM

Manufactured by:
Beijing Choice Electronic Technology Co., Ltd.
2nd and 3rd floors, No.2 Building, No.9 Shuangyuan Road,
Shijingshan District,100041 Beijing, PEOPLE'S REPUBLIC OF CHINA

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
IMPORTANT


Some users may experience delays while trying to obtain a reading. For some, numbers may show instantly, while for others, they may need to keep their finger in the oximeter longer. If you are not getting a reading after 30 seconds, please warm your fingers to increase circulation and try again.


The oximeter readings should not be used as the sole indicators of your wellness and exercise regimens. Use the readings in conjunction with other signs of excessive exercise, such as sore muscle pain and fatigue, to determine your daily training.

Quick Operation Guide

■ Before using the oximeter, please remove protective sticker covering the display.

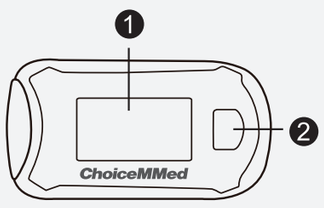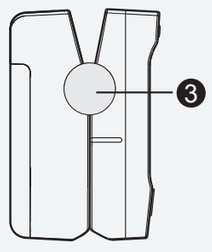
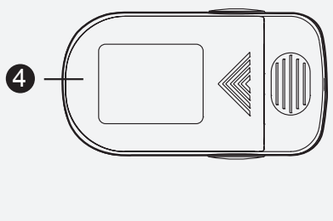
1. Displaying Screen
2. Power Key
3. Shaft Cover
4. Battery Box Cover

■ Brief Description of Front Panel

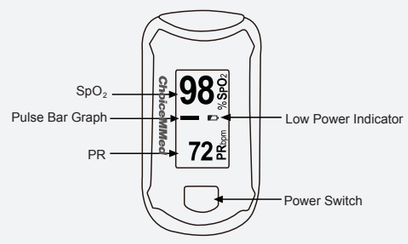The pulse bar graph displays corresponding with the patient's pulse beat. The height of the bar graph shows the patient's pulse strength.

■ Battery Installation

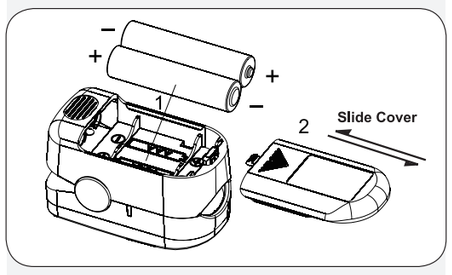
Please put in the batteries according to the polarity marked inside the battery compartment.

■ Using the Lanyard

1. Thread thinner end of the lanyard through the loop.
2. Thread thicker end of the lanyard through the threaded end before pulling it tightly.


Warnings!

- Keep the oximeter away from young children. Small items such as the battery door, battery, and lanyard are choking hazards.
- Do not hang the lanyard from the device's electrical wire.

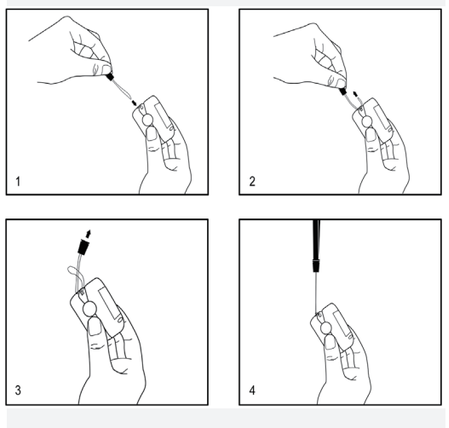
■ Operation Instruction

Place one of your fingers into the pulse oximeter to the end and press the switch button one time on front panel to turn it on.

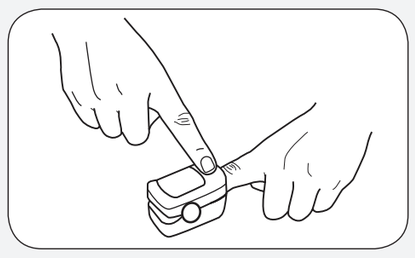There are 6 display modes. After turning on the pulse oximeter, each time you press the power switch, the pulse oximeter will switch to another display mode.

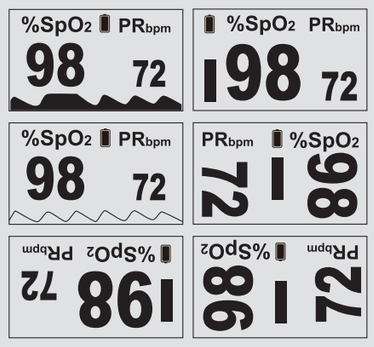When you press the power switch for more than one second, the brightness of the pulse oximeter will be changed by degrees, indicated by Br 1 and Br 10 in the display. There are 10 levels on brightness. The default level is level 4.

Keep your hands still for the reading.
The pluse oximeter will power off automatically in 8 seconds if there's no finger inside.

■ Warnings and Notes

Warnings:
1. Keep the pulse oximeter away from young children. Small parts such as the battery door and the batteries etc. May be hazardous if swallowed.
2. The lanyard may cause strangulation in conditions that may cause it to twist around the neck.

Notes:
1. Read the manual carefully before use.
2. The illustration used in this manual may differ slightly from the appearance of the actual product.
3. Follow local ordinances and recycling instructions regarding disposal or recycling of the device and device components, including batteries.

■ Frequently Asked Questions

What is a Pulse Oximeter?

A pulse oximeter is a non-invasive device that indirectly monitors blood oxygen saturation (SpO_2) and pulse rate (heart rate). It displays both blood oxygen saturation (SpO_2) and pulse rate (heart rate). Pulse oximeters provide an easy way of assessing your blood oxygen level and pulse rate.

What is SpO_2?

SpO_2 is also known as oxygen saturation. Oxygen saturation is a measure of how much oxygen the blood is carrying as a percentage of the maximum it could carry.

What is the normal range of SpO_2?

The normal range for SpO_2 is typically considered from 95%~99%. The SpO_2 measurement may be lower for people who live at high altitudes. Ask your health professional this question as it pertains to you.

What is the normal range for pulse rate?

The normal resting range for pulse rate is typically considered from 60~100 beats per minute. Ask your health professional this question as it pertains to you.

What is the waveform?

Each wave on the waveform is a heartbeat. During exercise and relaxation you can visually see your heartbeats increasing or decreasing respectively.

What kind of conditions may cause an inaccurate reading?

Cold hands, poor circulation, very weak pulse, movement, fingernail polish and acrylic nails may cause inaccurate results.

The SpO_2 is not changing – it's stuck?

SpO_2 does not change like pulse rate. It is slow to change.

The pulse rate is changing rapidly.

Your heart rate changes with emotions, excitement and exercise.

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© 2025 ChoiceMMed. ALL RIGHTS RESERVED