BIOS℠
Diagnostics

Cardiology Stethoscope 508DI



Read instructions carefully before use.

BIOS℠
Diagnostics

508DI Cardiology Stethoscope

Thank you for purchasing the BIOS Diagnostics™ Cardiology Stethoscope. Please read these instructions thoroughly and carefully before attempting to use this product.

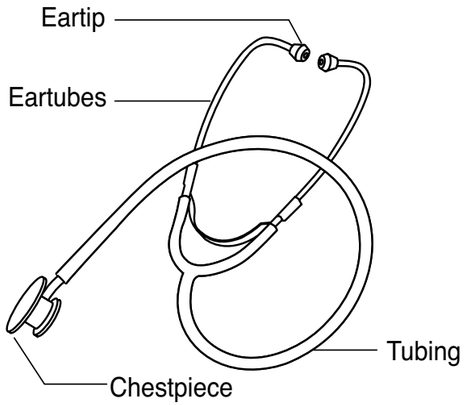


fig.1

This stethoscope has a dual head to hear high or low frequencies, simply rotate the chestpiece 180°.

⚠The stethoscope must be placed on the desired area of the body using only reasonable pressure.

Assemble and Use

1. Select the eartip that you prefer. For maximum acoustic seal and comfortable fit, the cardiology stethoscope is equipped with comfort sealing eartips. The comfort sealing tip is universal in size, and designed to conform to the ear canal so as to block out ambient noise. The eartips are screwed onto the tips for a secure connection.

Note: Different sizes of eartips are included. Simply remove eartips from stethoscope and attach the size that's more suitable.

2. To adjust the angle of the eartubes, hold the tubing in one hand and rotate each eartube to the correct angle as shown in fig.2.

3. To reduce spring tension, grasp each eartube at its bend and gently pull outwards as shown in fig.3.

4. To increase spring tension, squeeze the eartubes as shown in fig.4.

5. When the eartips are in your ears, the eartubes should be pointed forward as shown in fig.5. Improperly wearing the stethoscope headset can result in poor acoustical seal, and in some cases, complete sound blockage.

6. Only one side of the cardiology stethoscope works at one time. You need to open, or close the bell or membrane by rotating the chest piece. If the membrane is open, the bell will be closed, preventing sound from coming in through the bell, and vice versa.

⚠ Make sure that the eartips are pointed slightly forward into the ear canal.

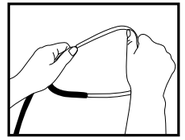


fig.2

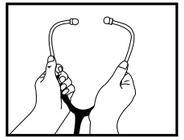


fig.3

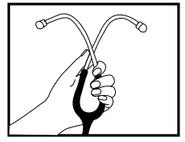


fig.4



fig.5

Cleaning

In order to avoid dirt accumulation and dust, always store the stethoscope in its packaging. The outside of the instrument can be cleaned with a damp, soft and fluff-free fabric. For disinfection, the stethoscope can be cleaned with fabric moistened with alcohol.

⚠ Do not use thermal sterilization or soaking.

Warning

Never bring the stethoscope into contact with liquid, and make sure that no liquid penetrates the housing! Since the stethoscope is not designed for operations, simple cleaning and/or disinfection with an alcohol based surface disinfectant is sufficient.
⚠ Do not use any scrubbing cleansing agents!

5 Year Warranty

The Cardiology Stethoscope has a 5 year warranty to be free of manufacturing defects in materials and workmanship under normal applications of the original owner. If this product becomes inoperable due to defect, return the product with all component pieces and proof of purchase to the address listed below. This warranty does not cover any shipping/transport costs. This warranty does not apply if the product is subject to misuse, neglect, rough handling or damage.

Ship the unit prepaid and insured (at owner's option) to:

Thermor Ltd.
Attn: Returns Department
16975 Leslie Street
Newmarket, ON L3Y 9A1
www.biosmedical.com
Email: support@biosmedical.com

BIOS Medical Hotline: 1-866-536-2289

FSC
www.fsc.org
FPO

Paper
Papier

BIOS | Medical
MANUFACTURED BY/FABRIQUÉ PAR :
THERMOR LTD.
16975 LESLIE STREET
NEWMARKET, ON L3Y 9A1
MADE IN CHINA/FABRIQUÉ EN CHINE
WWW.BIOSMEDICAL.COM
100% Canadian Owned & Operated
Propriétaire-exploitant canadien à 100 %
V. C