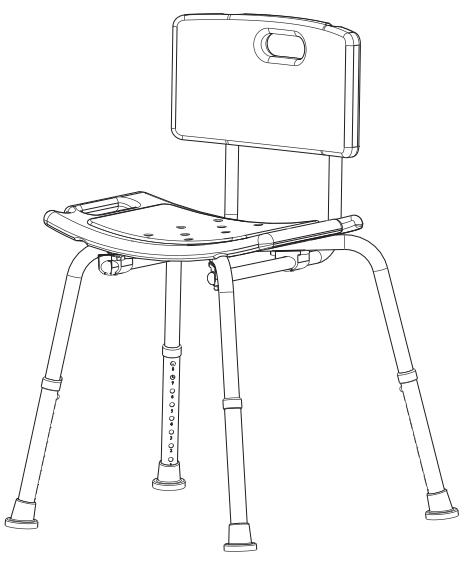
BIOS living™

59001
ADJUSTABLE BATH BENCH with backrest



Thank you for purchasing the BIOS Living™ Adjustable Bath Bench with Backrest. Please read all instructions carefully before attempting to install your Bath Bench.

CAUTION

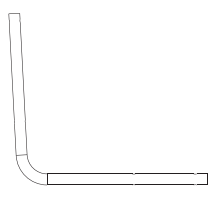
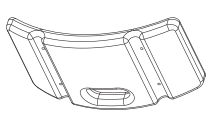
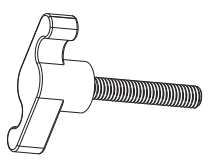
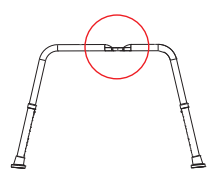
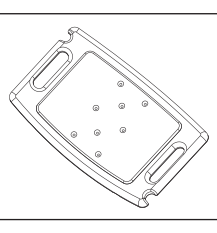
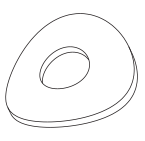
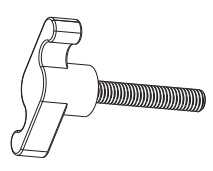
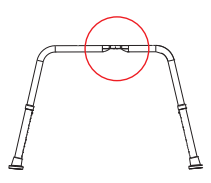
Ensure that the bench is correctly and securely locked in place before using.

Make sure to tighten the screws before each use and check that the height of all four legs is in the same numbered push pin receptacle.

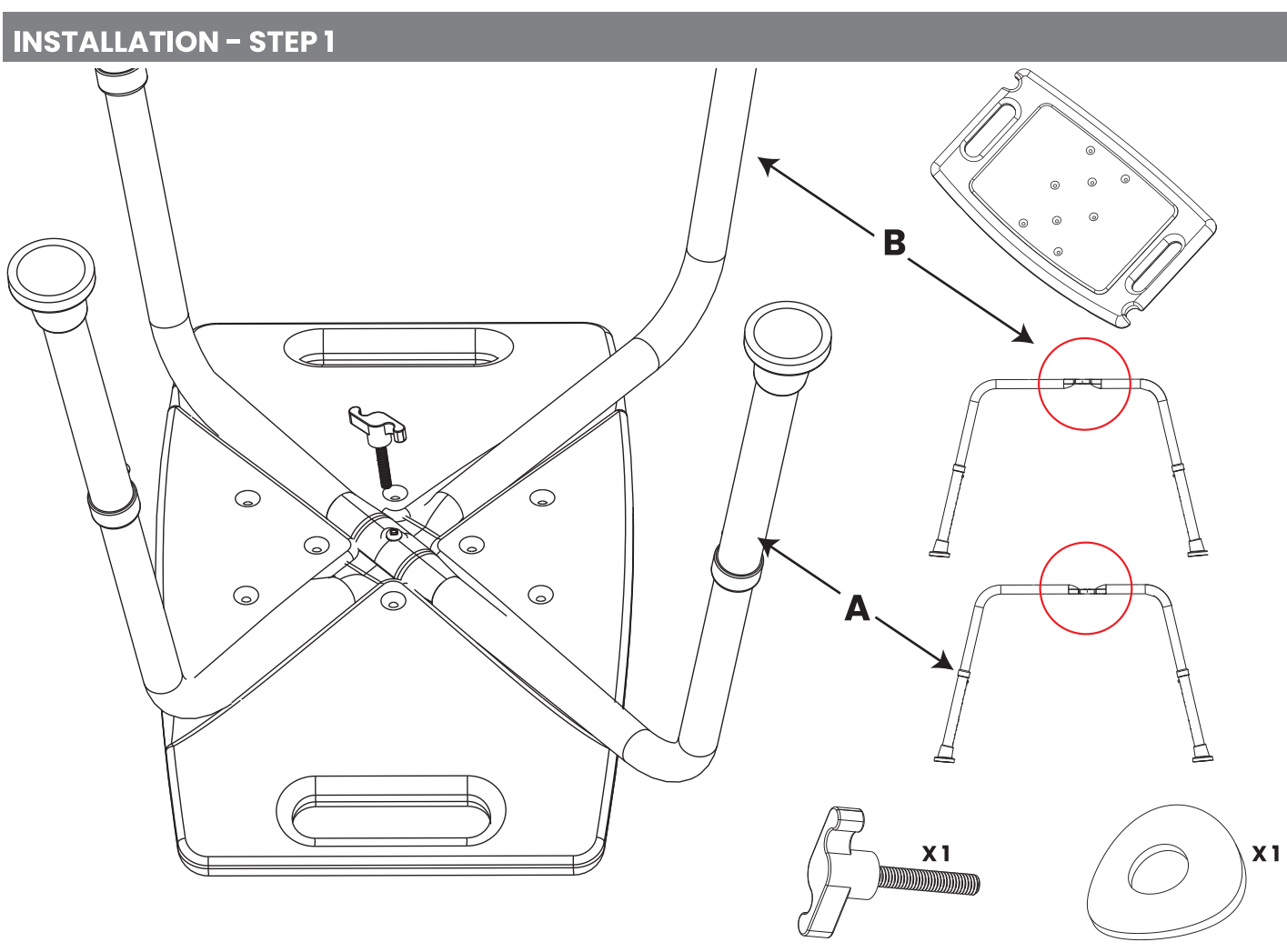
- **Do not** attempt to assemble the BIOS Living™ Bath Bench beyond your physical capability and/or without assistance.
- Consult with your healthcare provider to determine the suitability of the product for the particular user.
- **Do not** exceed the maximum weight limit of 300 lbs/136 kg. Doing so may result in injury.
- The BIOS Living™ Bath Bench should not be used as a step.

Note: the owner's manual was compiled from the most recent specifications and product information available at the time of the publication. Thermor Ltd. reserves the right to make changes, which may cause slight variations between the illustration and explanations in this manual and the product purchased.

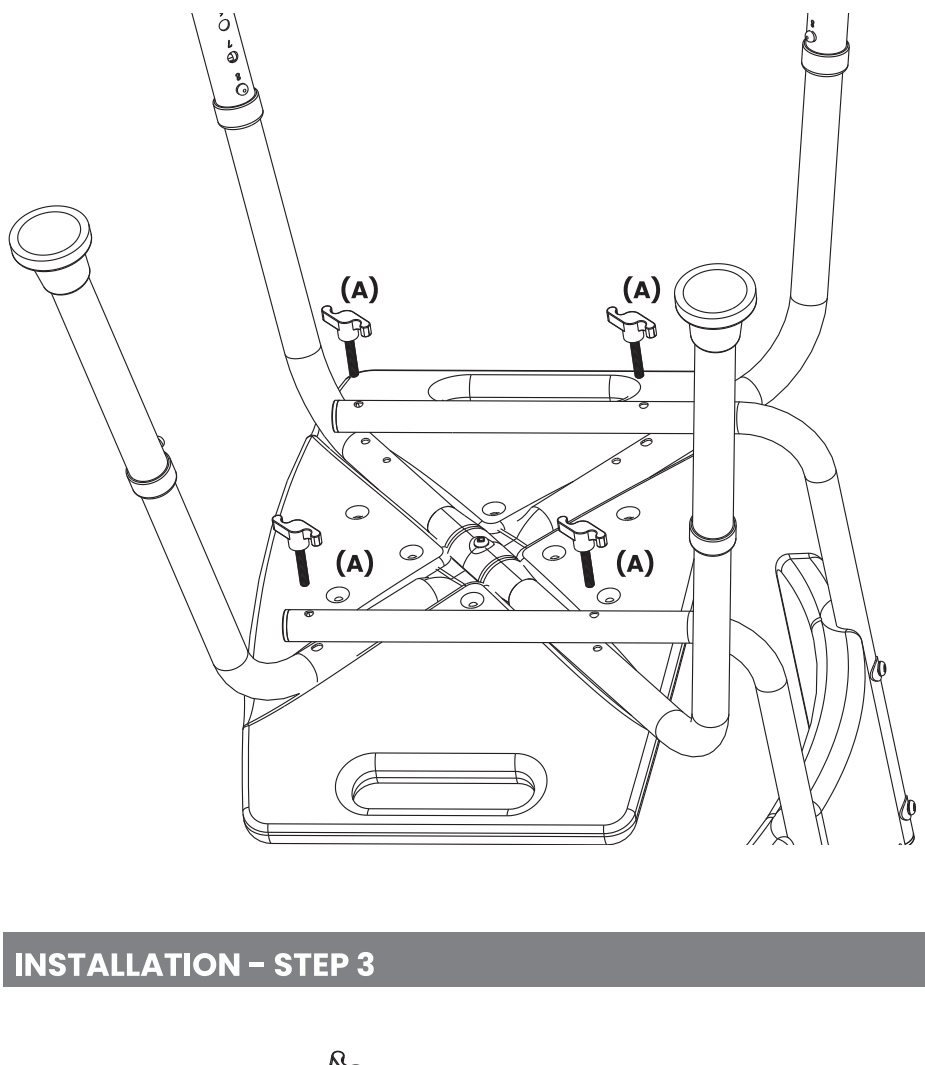
CHECK ALL PARTS

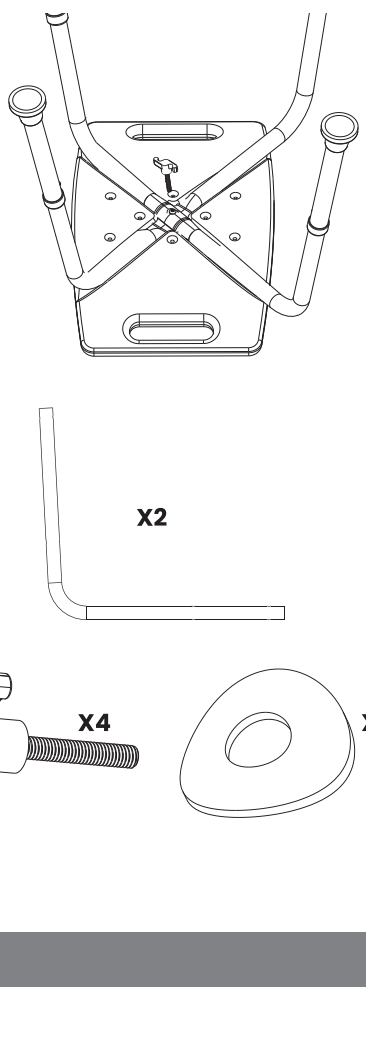
	x2 Backrest Tubes		x1 Plastic Backrest		x4 Knobs (55 mm) **(A)**
	x1 Frame Tube **(A)**		x1 Plastic Bench		x9 Arc Gaskets
	x5 Knobs (35mm) **(B)**		x1 Frame Tube **(B)**		

INSTALLATION – STEP 1

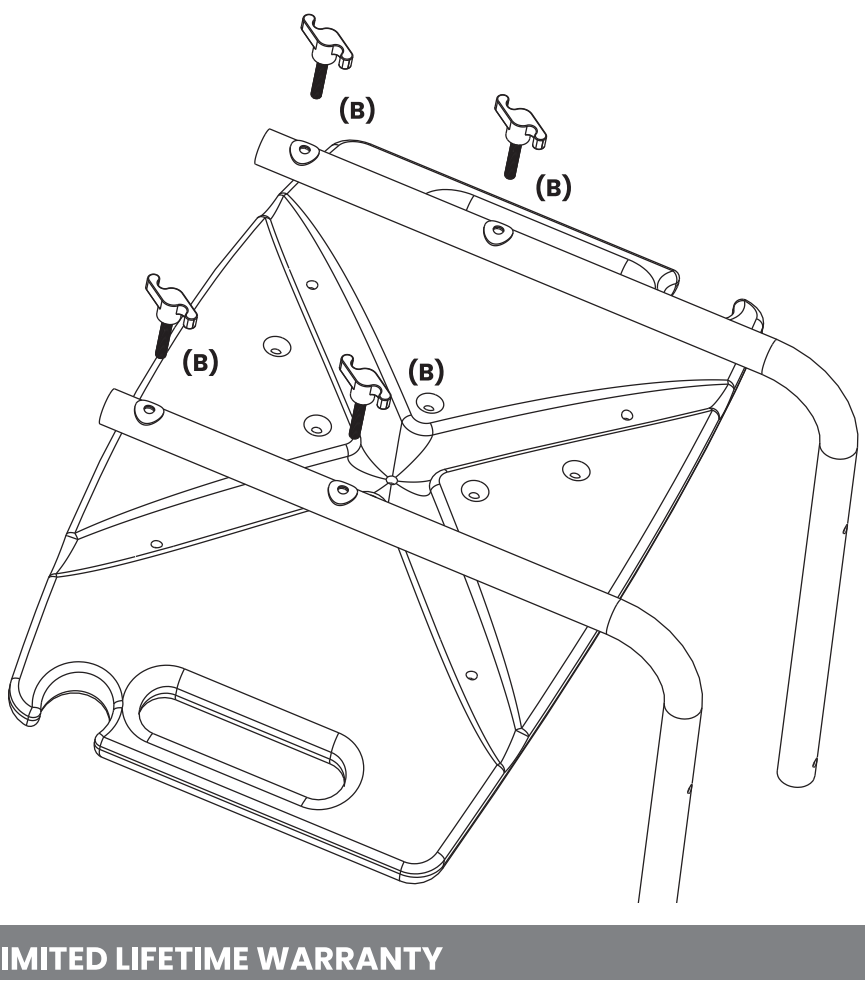


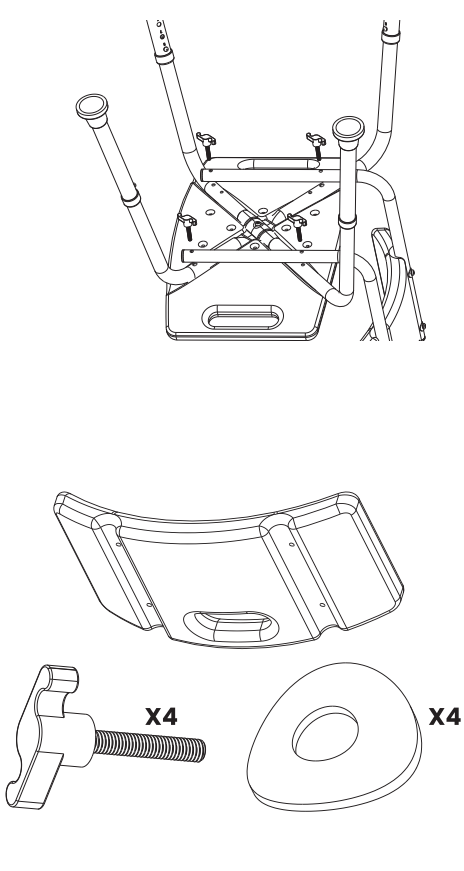
INSTALLATION – STEP 2





INSTALLATION – STEP 3





LIMITED LIFETIME WARRANTY

The BIOS Living™ Adjustable Bath Bench with Backrest has a Limited Lifetime Warranty to be free of manufacturing defects in materials and workmanship under normal applications, for the original owner. If this product becomes inoperable due to defect, return the product with all components and proof of purchase to the address listed below. This warranty does not cover any shipping/transport costs. This warranty does not apply to normal wear and tear or if the product is subject to misuse, neglect, rough handling or damage.

Ship the unit prepaid and insured (at owner's option) to:

Thermor Ltd.
Attn: Returns Department
16975 Leslie Street
Newmarket, ON L3Y 9A1
www.biosmedical.com
Email: support@biosmedical.com

FSC www.fsc.org FPO

Paper Papier

BIOS | Medical
MANUFACTURED BY/FABRIQUÉ PAR :
THERMOR LTD.
16975 LESLIE STREET
NEWMARKET, ON L3Y 9A1
MADE IN CHINA/FABRIQUÉ EN CHINE
WWW.BIOSMEDICAL.COM
100% Canadian Owned & Operated
Propriétaire-exploitant canadien à 100 %