BIOS living™

Adjustable Shower Stool
59005


Thank you for purchasing the BIOS Living™ Adjustable Shower Stool. Please read all instructions carefully before attempting to install your Shower Stool.

⚠ **CAUTION**


Ensure the seat is assembled correctly, bolts are tightened and stool placed on a flat surface securely before using

- **Do not** attempt to assemble the BIOS Living™ Adjustable Shower Stool beyond your physical capability and/or without assistance
- Consult with your healthcare provider to determine the suitability of the appliance for the particular user
- **Do not** exceed the maximum weight limit of 300 lbs / 136 kg – doing so may result in injury
- The BIOS Living™ Adjustable Shower Stool is not suitable for use as a step


STEP 1: CHECK ALL PARTS

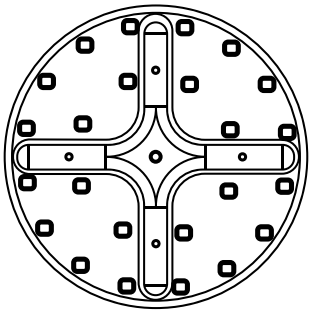


1 plastic seat

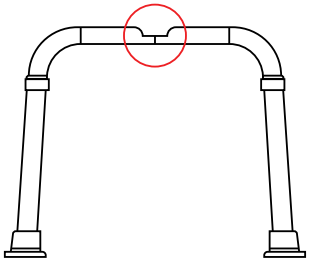


1 frame tube A

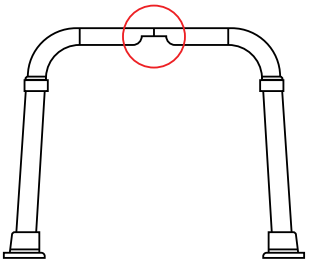

1 frame tube B

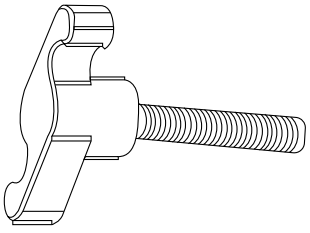


5 knobs

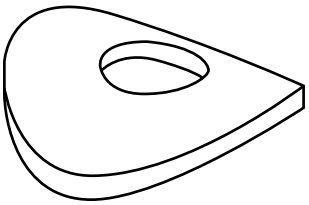


5 arc gaskets

STEP 2: ASSEMBLE

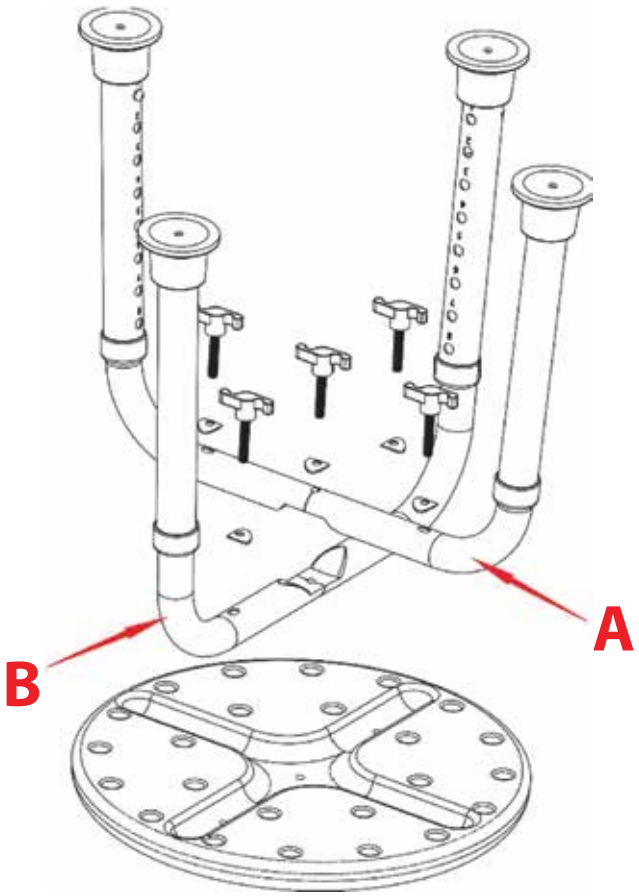

STEP 3: ADJUST FOR HEIGHT

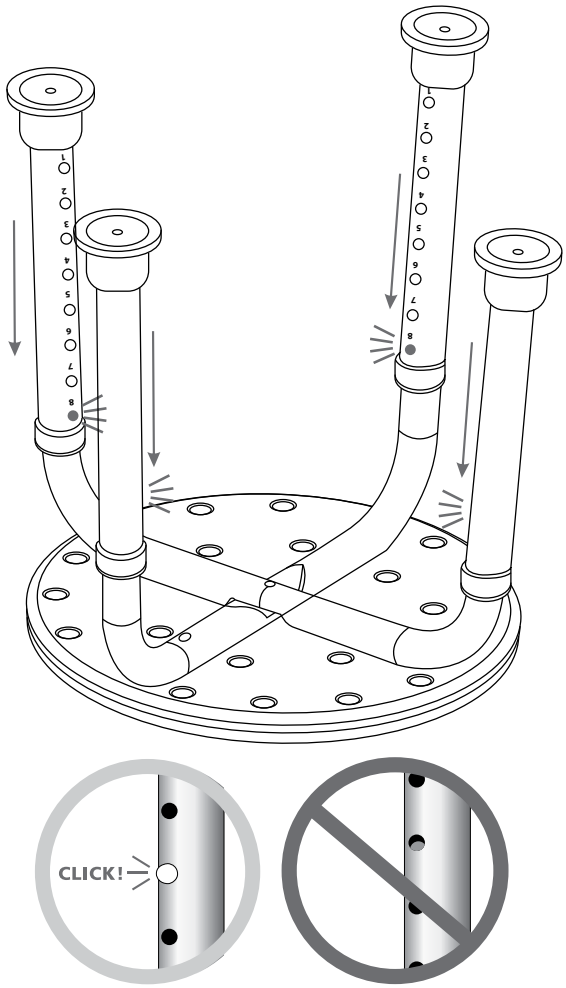

LIMITED LIFETIME WARRANTY

The BIOS Living™ Bath and Shower Stool has a Limited Lifetime Warranty to be free of manufacturing defects in materials and workmanship under normal applications, for the original owner. If this product becomes inoperable due to defect, return the product with all components and proof of purchase to the address listed below. This warranty does not cover any shipping/transport costs. This warranty does not apply to normal wear and tear or if the product is subject to misuse, neglect, rough handling or damage.

Ship the unit prepaid and insured (at owner's option) to:

Thermor Ltd.
Attn: Returns Department
16975 Leslie Street
Newmarket, ON L3Y 9A1
Email: support@biosmedical.com
Telephone: 1.866.517.2970

Paper
Papier

BIOS | Medical
MANUFACTURED BY/FABRIQUÉ PAR :
THERMOR LTD.
16975 LESLIE STREET
NEWMARKET, ON L3Y 9A1
MADE IN CHINA/FABRIQUÉ EN CHINE
WWW.BIOSMEDICAL.COM
100% Canadian Owned & Operated
Propriétaire-exploitant canadien à 100 %