
BIOS living™

Deluxe Commode 59024

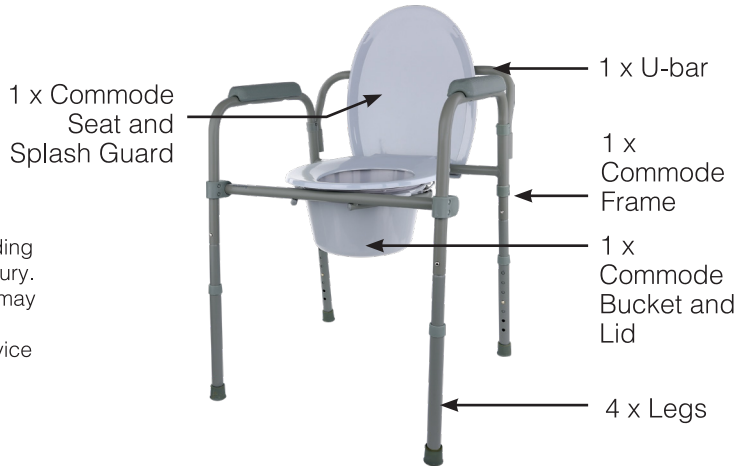

Caution

- **Do not** assemble or use this equipment without first reading and understanding this instruction manual. Failure to apply the cautions below could result in injury.
- **Do not** exceed the maximum weight capacity of 300 lbs / 136 kg. Doing so may result in injury.
- Consult with your healthcare provider in determining the suitability of the device for the particular user.

IMPORTANT: Always make sure all 4 legs are the same height before using.

Instructions



Unfold the commode frame.



Insert the front cross bar into the front brackets (located on the front legs of the commode frame). Ensure both metal pins "CLICK" completely into holes.





Adjust all 4 legs by moving the metal legs UP or DOWN over frame. Line up the metal pins with the corresponding holes; ensure the metal pins "CLICK" completely into holes. Make sure the height of each leg is the same.





"Snap" on the toilet seat and cover to the back support bar.



Attach the U-bar by inserting both ends into the back of the commode frame. Ensure the metal pins "CLICK" completely into holes.



Use the commode with the bucket for bedside and the splash guard with the commode for over the toilet use.

Folding the Commode

- Adjust legs to the lowest setting or remove to optimize space (optional)
- Push-in the metal pins and remove the U-bar
- Remove bucket or splash guard
- Push-in the metal pins to release the front cross bar and fold it down
- Fold down the toilet seat and lid
- Fold the commode as shown



Warranty

The Commode has a 3 year warranty on the frame and a 1 year warranty on plastic components (toilet seat, lid, bucket and splash guard) to be free of manufacturing defects in materials and workmanship under normal applications for the original owner. If this product becomes inoperable due to defect, return the product with all component pieces and proof of purchase to the address listed below. This warranty does not cover any shipping/transport costs. This warranty does not apply if the product is subject to misuse, neglect, rough handling or damage.

Ship the unit prepaid and insured (at owners option) to:
Thermor Ltd.
Attn: Returns Department
16975 Leslie Street
Newmarket, ON L3Y 9A1

www.biosmedical.com
Email: support@biosmedical.com

Tool Free Assembly

BIOS | Medical
MANUFACTURED BY/FABRIQUÉ PAR :
THERMOR LTD.
16975 LESLIE STREET
NEWMARKET, ON L3Y 9A1
MADE IN CHINA/FABRIQUÉ EN CHINE
WWW.BIOSMEDICAL.COM
100% Canadian Owned & Operated
Propriétaire-exploitant canadien à 100 %