BIOS living™

TOILET SAFETY FRAME INSTRUCTION MANUAL LF759

- Lightweight aluminum frame provides sturdy support
- Contoured armrests provide comfort and security
- Adjustable height: 28" - 32"/71 - 81 cm
- Adjustable width: 18" - 25"/46 - 64 cm
- Weight capacity: 300 lbs/136 kg

ASSEMBLY INSTRUCTIONS:

Step 1. Remove the toilet seat.
Begin by first removing toilet seat bolts.

Step 2. Align the mounting bracket with the toilet seat bolt holes. **(Figure 1)**

Step 3. Place the toilet seat on top of the mounting bracket and align holes.

Step 4. Install bolts and tighten securely.

Step 5. a) Slide armrest tubing into the mounting bracket and secure with push buttons. **(Figure 2)**

b) Assemble safety frame by inserting short tube into arm tube by sliding into arm rest and secure with "E" clip. (**Figure 2**)

Step 6. Set safety frame to desired size using push buttons. Make sure all buttons are secured before using the product. **(Figure 3)**

THREE YEAR LIMITED WARRANTY

The Toilet Safety Frame has a 3 year warranty to be free of manufacturing defects in materials and workmanship under normal applications for the original owner. If this product becomes inoperable due to defect return the product with all component pieces and proof of purchase to the address listed below. This warranty does not cover any shipping/transport costs. This warranty does not apply if the product is subject to misuse, neglect, rough handling or damage.

Thermor Ltd.
Returns Department
16975 Leslie Street
Newmarket, ON L3Y 9A1

www.biosmedical.com
Email: support@biosmedical.com
Telephone: 1.866.517.2970

BIOS | Medical
MANUFACTURED BY/FABRIQUÉ PAR :
THERMOR LTD.
16975 LESLIE STREET
NEWMARKET, ON L3Y 9A1
MADE IN CHINA/FABRIQUÉ EN CHINE
WWW.BIOSMEDICAL.COM
100% Canadian Owned & Operated
Propriétaire-exploitant canadien à 100 %
V. D

Paper

Papier

Figure 1

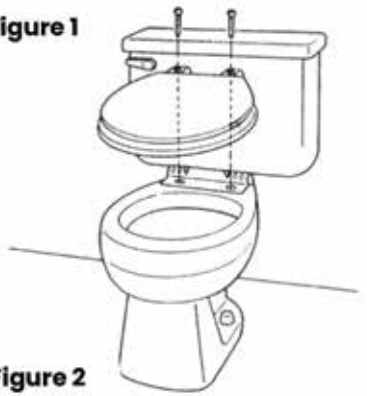


Figure 2

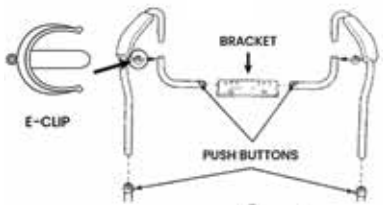


Figure 3

