BIOS living

LF775 & LF776
Hinged Toilet Riser Instructions

INTENDED USE


Raised toilet seats increase the independence of the user by adding height to most toilets. By easing the process of sitting and rising from the toilet, raised toilet seats provide comfort as well as support. The special feature of a hinged toilet riser is that it allows you to lift the riser and your toilet seat for cleaning.

NOTE: This product has been designed to fit most toilets. Ensure that the unit is correctly installed, fits properly and is stable before use.

CARE: We recommend using mild soap and water or anti bacterial wipes for cleaning.

Tools needed : flathead screwdriver and wrench

CONTENTS

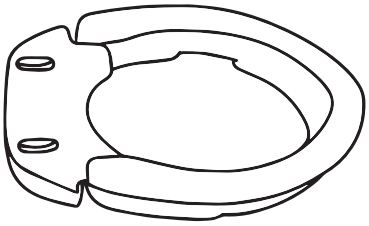

1 x Hinged Toilet Seat




2 x Bolts



4 x Metal Washers



4 x Plastic Washers

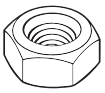

2 x nuts

(Steps 1-4, refer to Image A)

1. Remove your toilet seat from your toilet.
2. Position the hinged toilet riser on your toilet.
3. Place your toilet seat on top of the riser.
4. Line up the holes with the commode, slide the metal washer onto the bolt then thread through the toilet seat, riser and commode.

(Steps 5-8, refer to image B)

5. On the underside of the toilet use the plastic washer and screw on the nut by hand.
6. Tighten the bolt using a wrench to hold the nut and a flathead screwdriver to turn the bolt.
7. Tighten each side alternately until stable.
8. Check regularly that the toilet seat and riser are stable and secure before use. Tighten when needed.

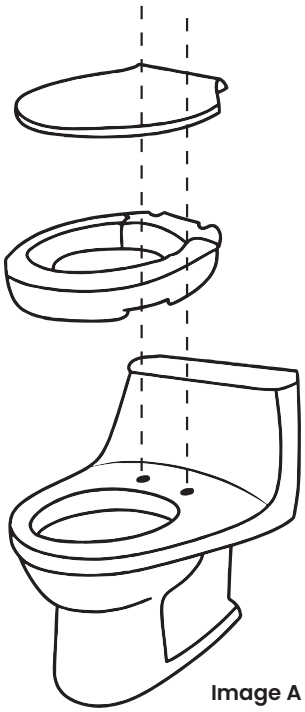


Image A

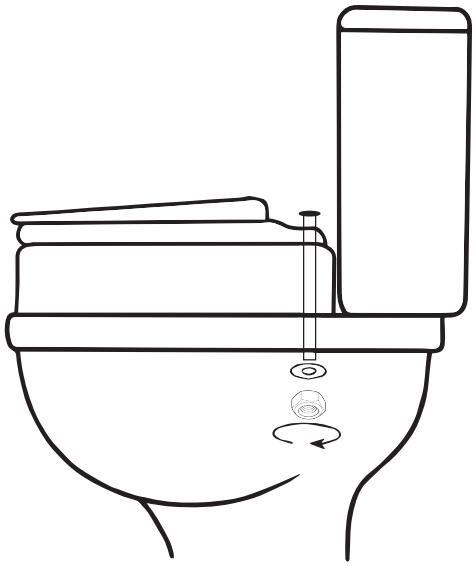


Image B

1 YEAR LIMITED WARRANTY

The Bios Living Hinged Toilet Seat has a limited 1 year warranty to be free of manufacturing defects in materials and workmanship under normal applications, for the original owner. If this product becomes inoperable due to defect, return the product with all components and proof of purchase to the address listed below. This warranty does not cover any shipping/transport costs. This warranty does not apply to normal wear and tear or if the product is subject to misuse, neglect, rough handling or damage.

Ship the unit prepaid and insured (at owner's option) to:

Thermor Ltd.
Attn: Returns Department
16975 Leslie Street
Newmarket, ON L3Y 9A1

Paper
Papier

FSC
www.fsc.org
FPO

BIOS | Medical
MANUFACTURED BY/FABRIQUÉ PAR :
THERMOR LTD.
16975 LESLIE STREET
NEWMARKET, ON L3Y 9A1
MADE IN CHINA/FABRIQUÉ EN CHINE
WWW.BIOSMEDICAL.COM